

# SRG-3

## Spinning Rotor Gauge System



The Spinning Rotor Gauge (SRG) line is well established with users around the world as a stable, high accuracy transfer standard for high vacuum pressure measurements. We continue to be at the forefront of SRG development with the introduction of the SRG-3 Spinning Rotor Gauge System.

The complete SRG-3 system is comprised of a microprocessor-based electronics unit, a measuring head and a ball/flange assembly. The measuring head easily decouples from the ball/flange assembly, allowing several permanently mounted ball/flange assemblies to be used with one measuring head and electronics unit. This facilitates taking measurements from multiple locations when continuous readings are not necessary.

### Product Features

- German Calibration Service, DKD traceable
- Recognized transfer standard
- Guaranteed precision from  $5 \times 10^{-7}$  to  $1 \times 10^{-2}$  mbar (Torr)
- Overall measuring range extends to 1 mbar (Torr)
- No thermal or ionization effects
- Wetted surfaces are stainless steel



### Key Benefits

- Ensures absolute consistency
- Superior long term stability
- Assures chamber-to-chamber consistency

## Principle of a Spinning Rotor Gauge


The sensing element of the SRG is a magnetically levitated stainless steel ball that resides in a tube connected to the vacuum system. The ball is suspended and rotated by the separable drive assembly. Once the ball reaches a certain rotational speed, the ball is allowed to free spin. Collisions of the gas molecules with the surface of the ball impart a drag which decelerates the ball. Pickup coils measure the deceleration. The electronics unit then uses an algorithm based on the sampling time, physical parameters and gas characteristics to accurately determine the pressure.

Since the measurement principle is mechanical, there are no thermal or ionization effects (as is the case with ionization gauges) that will influence the accuracy of the measurements. Also, the calibration constant is only weakly influenced by the ball's surface roughness, making the SRG a very stable reference standard.

The SRG-3 will find applications in metrology and in experiments where high accuracy and repeatability are essential. Away from the test bench, the MKS Spinning Rotor Gauge may also be used in situ on a process chamber to assure high vacuum gauge reading accuracy and tool-to-tool consistency. This provides repeatable process results, high uptime and optimum cycle times.


MKS Instruments, Inc. delivers the SRG-3 with a calibrated measuring head and ball/flange assembly along with a calibration certificate. Calibration is made against the reference standards of our DKD1 calibration laboratory. As a result, the Spinning Rotor Gauge SRG-3 system offers the user consistent traceability as required for ISO9000.

Should measuring flanges with non-calibrated balls be ordered, they will be checked for usability by the factory.

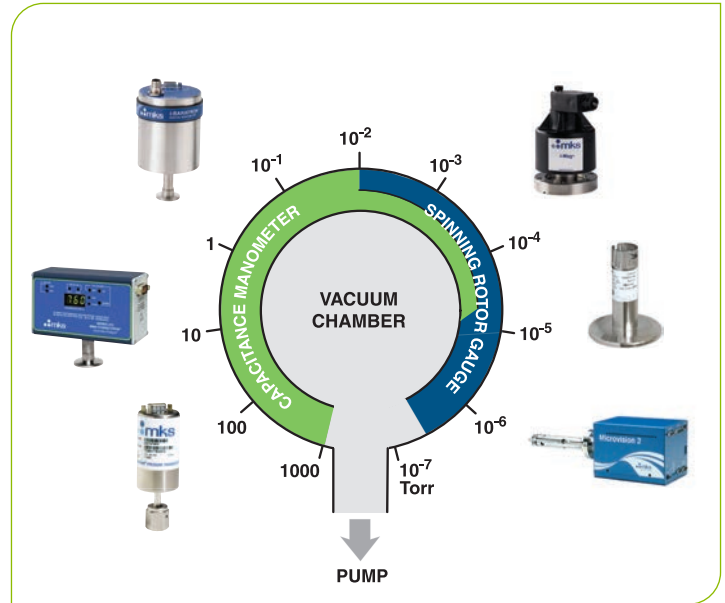

### Calibration, Traceability, Stability

In addition to the SRG-3, MKS Instruments offers transfer standard Baratron® capacitance manometers to meet all of your needs for vacuum gauge calibration.

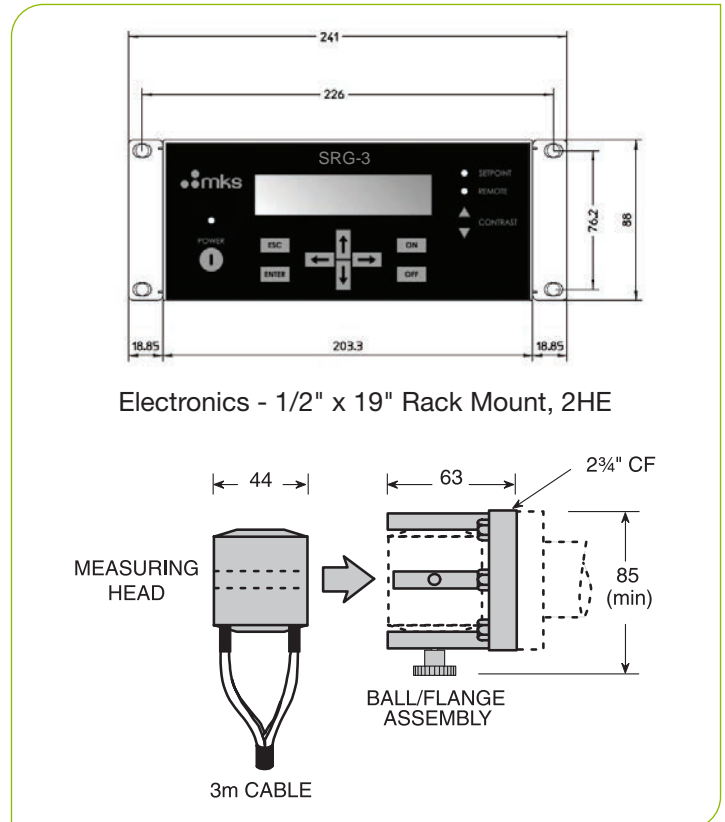

### Dimensional Drawing

Note: Unless otherwise specified, dimensions are nominal values in millimeters.

# Specifications

<b>Measuring Range</b>	5x10 <sup>-5</sup> to 100 Pa
<b>Accuracy</b>	
Up to 1 Pa	1% of measuring value + U; U expresses the residual drag variation
1 to 100 Pa	Increasing up to 10% of measured value (typical)
<b>Long-Term Stability</b>	Better than 1% per year
<b>Display</b>	4½ digit LCD display
<b>Pressure Display</b>	Pa, mbar, or Torr; user selectable
<b>Display Format</b>	Scientific notation
<b>Analog Output</b>	0-10 VDC, into 2kΩ min. load, linear or logarithmic
<b>Switching Points</b>	2 SPDT relays, adjustable, nominal switching capacity (resistive load): 1A 30VDC, 0.5A 25VAC eff
<b>Digital Interface</b>	(a) RS-232; (b) USB (USB 2.0 compatible)
<b>Measuring Flange</b>	2¾ CF, made of stainless steel, bakeable up to 450°C (without measuring head)
<b>Rotor</b>	Stainless steel ball, 4.5 mm in diameter (nominal)
<b>Operating Temperature</b>	
Electronics	10-40°C
Measuring Head	10-50°C
Ball/Flange Assembly	Bakeable to 450°C
<b>Enclosure</b>	Half 19-inch rack mount or table top; 241mm x 185 mm x 88 mm (BxDxH)
<b>Weight</b>	
Electronics	2.5kg
Measuring Head	0.6kg
Measuring Flange	0.5kg
<b>Operating Voltage</b>	85-256 VAC, 47-63 Hz
<b>Wetted Materials</b>	Ball and Flange
<b>Programming Features</b>	15 programmable data sets to minimize set up time; database of 17 frequently used gases with the option for the user to add 8 additional gases
<b>Compliance</b>	CE

# Ordering Information

---

## SRG-3 Spinning Rotor Gauge

---

SRG3-EL  
SRG-SH700-V3  
SRG-BF-CAL  
SRG-BF

Electronic Control Unit  
Measuring head with 3m long cable  
Ball/Flange Assembly, calibrated  
Ball/Flange Assembly, not calibrated

---



[www.MKSINST.com](http://www.MKSINST.com)

+1-978-645-5500 | +1-800-227-8766

SRG3\_05/20  
©2020 MKS Instruments, Inc.  
Specifications are subject to change without notice.

---

MKS products provided subject to the US Export Regulations. Export, re-export, diversion or transfer contrary to US law (and local country law) is prohibited.

mksinst<sup>®</sup> is a trademark and Baratron<sup>®</sup> is a registered trademark of MKS Instruments, Inc., Andover, MA.